



iCode[®]

QUICK REFERENCE SHEET

www.icodeconnect.com

www.bmc-iCode.com

What is iCode[®]?

iCode[®] allows remote access of sleep compliance data with a phone call.

- iCode[®] is web-based. No software installation or direct cable connections are necessary.
- iCode[®] is free, and comes standard on every model.
- Patient auto-dial with voice recognition for automated compliance data retrieval.



How does it work?

Sleep compliance data is encoded into a 16 character numeric string, which is displayed on the device screen. Using our online report generator at www.bmc-iCode.com or icodeconnect.com, the iCode[®] string is decoded. An easy-to-read report may be printed for placement in the record or saved as a pdf file.

Just a Phone Call Away

iCode[®] allows quick access to both compliance data needed for reimbursement and summary sleep data for therapy monitoring with just a phone call.


iCode[®] data is objective as required by the Medicare CPAP Local Coverage Determination (LCD). To prevent fraudulent use of a compliant patient's iCode[®] being used for another patient, both the iCode[®] string and the device serial

number are check digit validated. Entering an iCode[®] string from one CPAP device with the serial number of another CPAP device will result in an error. Moreover, iCode[®] is HIPAA compliant. No patient information is stored on the server after the report is generated.

Best 30 Day Adherence Score

iCode[®] helps with reimbursement by providing a Best 30 Day Adherence Score. The best 30 day is a percentage of the best consecutive 30 days of usage that meets the CPAP LCD. Using an adherence threshold of 4 hours or more of continuous use each night, the Best 30 Day Adherence Score would show the best adherence percentage achieved for any rolling period of 30 consecutive days of use within a 90 day time frame.

Interpreting Your iCode®

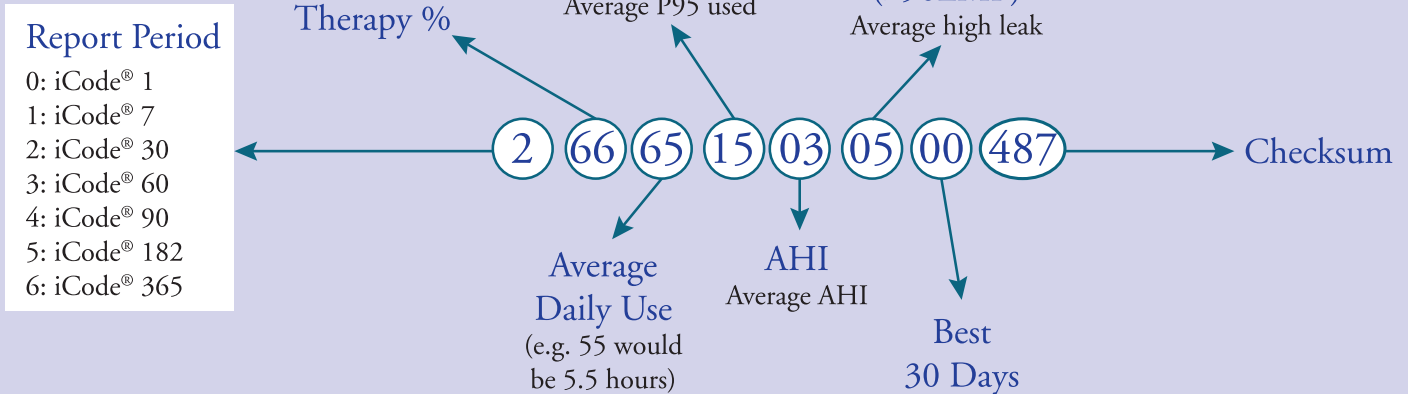

iCode Connect

<https://icodeconnect.com>

iCode Compliance Report for Smith, John

Date Collected: 3/04/2014

Patient information		Provider information	
Patient name:	John Smith	Organization:	Sample Co., Ltd.
Gender:	Male	Street 1:	1234 Maple St.
Date of birth:	4/04/1960	Street 2:	
Therapy start date:	3/04/2014	City:	Lake Wales
Device:	CTX1309J112	State:	FL
Patient ID:	JSmith1960	ZIP:	33853
Street 1:		Phone:	(555) 555-5555
Street 2:		Physician information	
City:	Lake Wales	Name:	
State:	FL	Organization:	
ZIP:	33859	Street 1:	
Report information		Street 2:	
Device serial number:	CTX1309J112	City:	
iCode string:	1716014020100668	State:	
		ZIP:	
		Phone:	





Generate a Sample Report

iCode® is easy to use. Generate a sample report using the steps above and the following sample iCode® strings:

Device Number: 1309J112

iCode® 1 : 0995514010100977
 iCode® 7 Days: 1716014020100668
 iCode® 30 Days: 2606514040300554
 iCode® 60 Days: 3676915050524249
 iCode® 90 Days: 4568017060624237
 iCode® 182 Days: 5777518080700049
 iCode® 365 Days: 6797220091000187

Report Date	iCode period	Best 30 days compliance	Days of therapy (>4h)	Avg. daily use time	Avg. P95	Avg. Pressure	Avg. AHI	Avg. SNI	High leak time (%)
3/04/2014	7	N/A	5 / 7 (71%)	6:00	7.0		2		1%
3/04/2014	90	24 / 30 (80%)	50 / 90 (56%)	8:00	8.5		6		6%
3/04/2014	182	N/A	140 / 182 (77%)	7:30	9.0		8		7%
3/04/2014	365	N/A	288 / 365 (79%)	7:12	10.0		9		10%

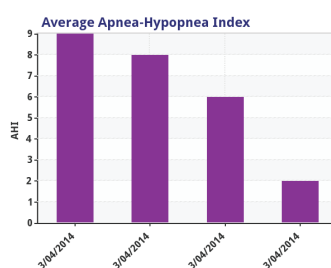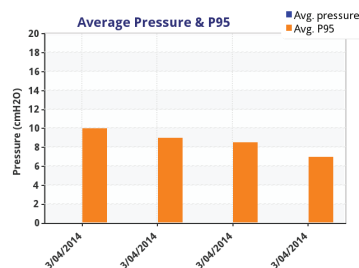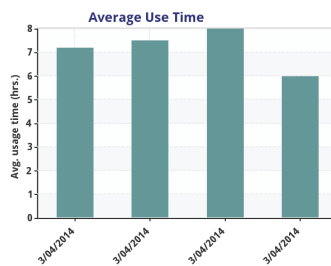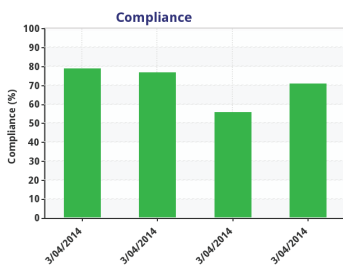

Generate and Archive iCode® Reports on iCode Connect™

For detailed instructions, see the user manual at www.icodeconnect.com.

iCode Connect™ allows you to easily and efficiently manage your patients' sleep therapy. Generate and store sleep reports that you can access any time and anywhere you have an internet connection. Set up a provider account today to get started.



www.3Bproducts.com

*Device serial numbers should be recorded and available in the patient record during set-up. Alternatively, the device serial number may be found on the metal plate on the bottom of the device.